



# AS608 / AS608e Screening Audiometers

---

Portable screening  
audiometer



  
**Interacoustics**

# Portable and easy to use

The AS608 / AS608e portable audiometer is ideal for schools, industry and primary medical practices or wherever quick hearing screening evaluations are performed. Easy to use, it requires minimal training prior to use and offers calibrated pure tone and warble tone stimuli over a wide range of frequencies.

## Portable

The AS608 / AS608e includes a dedicated, lightweight carrying bag with shoulder strap that will accommodate the AS608 / AS608e, headset, and audiogram charts.

## Batteries or power supply

The AS608 / AS608e will operate with 3 standard AA batteries or a medically approved external power supply. Typical battery life is six months.



The soft carrying case



Diagnostic Suite software (AS608e only)

## PC-integration (AS608e only)


The AS608e offers PC-integration through a USB cable using the Diagnostic Suite Windows® software for transferring audiometric data to the PC. The diagnostic suite offers features such as:

- Print Wizard
- Advanced report page
- Session handling
- Database storage: OtoAccess™ and Noah compatible


## Auto Threshold (AS608e only)

The AS608e includes the Hughson Westlake automatic pure tone test procedure. When the test is completed the results are easily recalled from the internal memory of the AS608e and displayed in the Diagnostic Suite PC software.



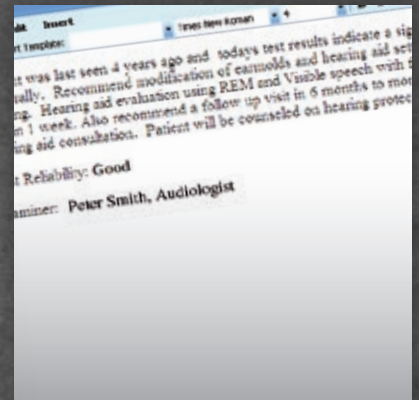

Report page from the Diagnostic Suite  
(AS608e only)

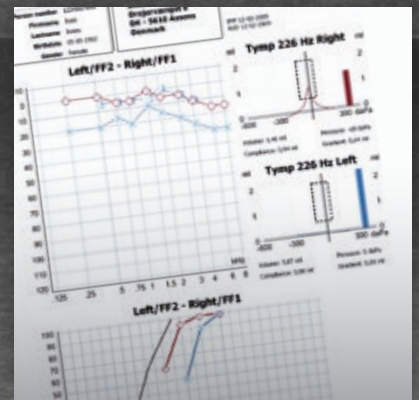

Print Layout from the Diagnostic Suite  
(AS608e only)

## A preview of the benefits

- Robust construction
- 11 frequencies from 125Hz to 8kHz
- Customize frequency test selections
- Output levels from -10dB to 100dBHL
- Pure tone, warble tone and pulsed tone
- Silent stimulus switch
- Battery status indicator
- Carrying bag
- Patient response switch (included with AS608e)
- External power supply (optional)
- Noise reducing headset (optional)
- Talk Forward (only AS608e)
- PC-integration (only AS608e)
- Hughson Westlake automatic pure tone test procedure (only AS608e)

**Total weight  
less than  
3.5 lbs**





---

## Interacoustics

Interacoustics is a world leading provider of diagnostic solutions in the field of hearing and balance assessment. Since 1967 we have designed and manufactured our innovative diagnostic solutions for the world of audiology with a constant focus on providing our customers with quality, dependable products.

---

## With you at all times

We operate in over 100 countries worldwide through a carefully selected network of distributors and service centres. Purchasing an Interacoustics product guarantees not only a quality product, but also direct access to our highly professional training and support service.

---

## Product specifications

All technical and hardware specifications concerning the AS608 can be downloaded from our website.

**Read more at**  
[interacoustics-us.com](http://interacoustics-us.com)

---

### Interacoustics USA

10393 West 70th Street  
Eden Prairie, MN 55344

**T** +1 800 947 6334

**F** +1 952 903 4200

[info@interacoustics-us.com](mailto:info@interacoustics-us.com)  
[interacoustics-us.com](http://interacoustics-us.com)

Go online to  
explore our  
full product  
range



# Interacoustics