

---

# Callisto™ Portable audiometry & REM

---

Redefining portability  
- Audiometry, REM, and  
Visible Speech Mapping  
in one box



# Fitting into a modern practice

Small footprint and lightweight (1.25 lbs)



The Callisto™ is a portable PC-based modular system that can include both a diagnostic audiometer and a full featured Real-Ear Measurement module (including Visible Speech Mapping).

## Highlights

- Intuitive - combines advanced testing and counseling
- Conforms to other Interacoustics Suite Software reducing the learning curve
- Mix and match Audiometry, REM and VSP
- Tailored customizable test routines
- Merge Video Otoscopy into the suite
- NOAH/OtoAccess™ compatible

## Bring the clinic to your client

- Ideal for travelling and in-home visits
- Small footprint and lightweight (1.25 lbs)
- Powered via the PC USB
- Optional dedicated carrying bag on wheels
- Leave Callisto™ in the bag while testing

## Real Ear to Coupler Difference

- Save fitting time, RECD feature will measure performance before the patient arrives
- The adult targets are not always compatible for children's smaller ear canal because greater SPL will be generated. Consequently, if using RECD when fitting a child there will be a lesser risk of over amplifying.




Visible Speech Mapping

## Audiometry

- Customizable test setups, printouts and report templates
- Visual and auditory counseling tools (Speech Banana, Hearing Loss Simulator, Master Hearing Aid)
- Digitized word lists from your PC hard-drive
- Audiogram historical comparison overlay
- Masking Help
- TEN test
- High Frequency capability and Digitized QuickSIN™ (optional)
- Customizable counseling overlays

## Real Ear Measure/Visible Speech Mapping

- Visible Speech Mapping for combined verification and counseling
- Multiple speech stimuli in several languages and live voice
- True 'On-Top' mode for fast interactive HA fittings
- NAL-NL2 and DSLv5
- Open Fit capable
- Directionality testing
- See data in graph or table mode
- Quickly switch between Response or Gain views
- In-the-ear monitor
- Binaural REM
- Compatible with AutoFit, REMfit and EXPRESSfit
- Genie module launch

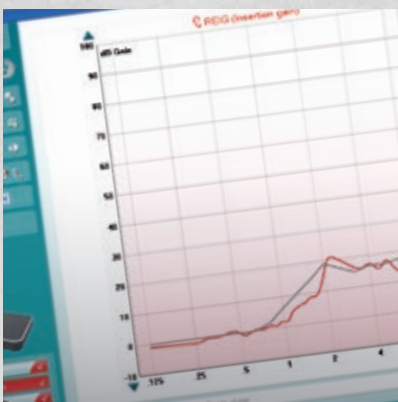

Real Ear Measurement (REM)

- On-top mode
- Open fit capability
- RECD and coupler verification
- DSLv5 and NAL-NL2 prescription targets

“Inspired by advanced Audiometry and Real-ear systems, but wrapped in a small package”

**Real Ear Measurements**  
- an Introduction:



**Visible Speech Mapping**  
- an Introduction:





## Interacoustics

Interacoustics is a world leading provider of diagnostic solutions in the field of hearing and balance assessment. Since 1967 we have designed and manufactured our innovative diagnostic solutions for the world of audiology with a constant focus on providing our customers with quality, dependable products.

## With you at all times

We operate in over 100 countries worldwide through a carefully selected network of distributors and service centres. Purchasing an Interacoustics product guarantees not only a quality product, but also direct access to our highly professional training and support service.

## Product specifications

All technical and hardware specifications concerning the Callisto™ can be downloaded from our website.

Read more at  
[interacoustics-us.com](http://interacoustics-us.com)

Go online to  
explore our  
full product  
range

### Interacoustics USA

10393 West 70th Street  
Eden Prairie, MN 55344

T +1 800 947 6334

F +1 952 903 4200

[info@interacoustics-us.com](mailto:info@interacoustics-us.com)  
[interacoustics-us.com](http://interacoustics-us.com)



# Interacoustics