

---

# The all new AA222


## A unique and powerful diagnostic combo

The AA222 combines  
the power of a clinical  
middle ear analyzer  
with a complete  
diagnostic audiometer.



**Interacoustics**


# A unique and powerful diagnostic combo

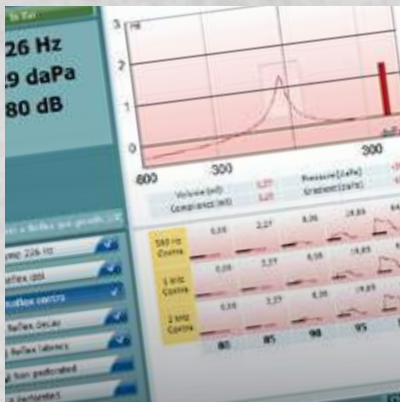

Tympanogram and reflex test results displayed simultaneously



Speech test available in both table and graphic mode

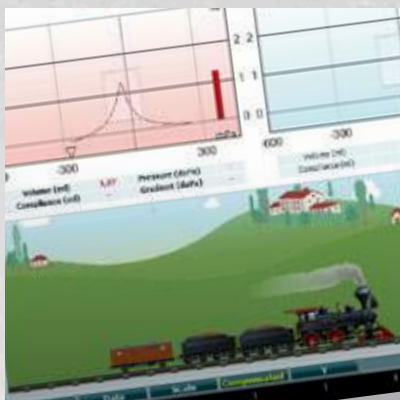

Child distractor - moving choo-choo train

The perfect combination device when space and cost savings need to be considered. The power of a clinical middle ear analyzer with a complete diagnostic audiometer makes the AA222 perfect for a start-up practice or satellite locations.

## Comprehensive test selections

- 226 Hz tympanometry (High Hz optional)
- Ipsi/contra reflexes
- Reflex decay
- ETF tests
- Pure tone air & bone conduction
- Speech testing
- Auto HL, Stenger, ABLB, Weber

## Outstanding features - impedance

- Child distractor - moving choo-choo train
- Large, 10" color display for easy viewing
- Durable clinical probe
- Optional diagnostic probe for fast screening applications
- Brilliant probe status lights
- Customized test sequences



## Outstanding features - audiometry

- Multiple user test setups
- Large, 10" color audiogram display
- Comprehensive library of speech recordings
- Binaural speech testing
- Special QuicksIN test mode
- Operator mic/monitor headset - or -
- Built-in gooseneck mic and speaker

## Outstanding performance

- Extremely fast tymp and reflex tests
- Customized setups for impedance and audiometry
- Direct printing to external printer
- Included database or NOAH interface
- Diagnostic Suite - merges test data into custom reports

# Tympanometry & audiometry in one powerful device



Color:	Probe:	Status:
Red		Right ear is selected. Probe is out of the ear.
Blue		Left ear is selected. Probe is out of the ear.
Green		Probe is in the ear and a seal is maintained.
Yellow		Probe is in the ear and is blocked or leaking.
White		The probe has just been attached. Probe status is yet unknown.

Probe status displayed on diagnostic probe / shoulder box for clinical probe



## Interacoustics

Interacoustics is a world leading provider of diagnostic solutions in the field of hearing and balance assessment. Since 1967 we have designed and manufactured our innovative diagnostic solutions for the world of audiology with a constant focus on providing our customers with quality, dependable products.

---

## With you at all times

We operate in over 100 countries worldwide through a carefully selected network of distributors and service centres. Purchasing an Interacoustics product guarantees not only a quality product, but also direct access to our highly professional training and support service.

---

## Product specifications

All technical and hardware specifications concerning the AA222 can be downloaded from our website.

Read more at  
[interacoustics-us.com](http://interacoustics-us.com)

---

### Interacoustics USA

10393 West 70th Street  
Eden Prairie, MN 55344

T +1 800 947 6334

F +1 952 903 4200

[info@interacoustics-us.com](mailto:info@interacoustics-us.com)  
[interacoustics-us.com](http://interacoustics-us.com)

Go online to  
explore our  
full product  
range



# Interacoustics

ABR/OAE

Balance Assessment

Fitting Systems

Middle Ear Analyzers

Audiometers